


# **KEN BWD 736 & BWD 738**

# Large capacity machines for hospitals



For more information consult – www.ken.dk

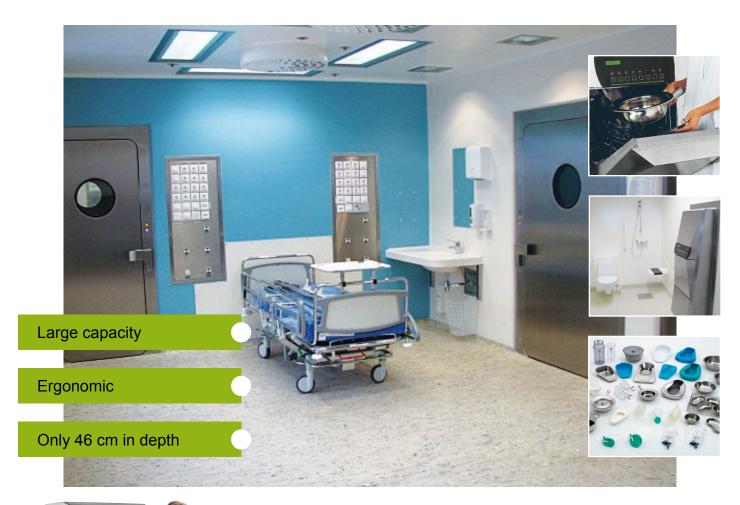

# Unique advantages with large capacity in minimum space

KEN BWD736 is available as wall or floor fitted, while BWD738 is a pass-through machine. Even with the door open the BWD is only 46cm in depth. Despite this the machine will hold 3 round bedpans with lid or 3 urinals.

The bedpan washer can be delivered with manual or automatic door operation.

By automatic door operation it is possible to select between a foot or remote solution. The remote control is placed on a wall so it can be activated by a slight pressure e.g. with the elbow by which means the loading and emptying are more efficient.

A wide variety of extra equipment ensures flexibility and optimum use.

www.ken.dk



### Advantages:

- Clear text screen displays information or error messages
- User-friendly operation panel with Proces Indication Bar
- Available with up to 5 programmes as standard
- Large wash chamber
- Possibility for SPRI basket 450-x340x70mm
- Up to 3 round bedpansor 3 urinals
- Pass-through model for isolation rooms

#### Service

All KEN machines conform to present EU Directives, and all hospital machines have been constructed in accordance with EN ISO 15883, including, for example, the testing of its cleaning capacity according to SIS-TR 3:2001 (formerly SPRI).

Contact your local KEN distributor to make sure you have the correct service for your machines. The right service will prolong the machine's working life.

#### **Accessories**

Look for equipment with the aid of the equipment selector at www.ken.dk, or contact KEN for information about further possible equipment, including customized special solutions.

### Contact a KEN specialist

KEN will gladly help you find the right solution. Consult www.ken.dk to find your KEN dealer.

Tlf.: +45 62 63 10 91 e-mail: ken@ken.dk
Internet: www.ken.dk



## **Specifications**

Power connection (standard)	3N~400V, N, J, 50Hz or multiple choices
Connected load/initial fuse	6kW/10A alternatively 9kW/16A
Water connection	Cold water: 1 pcs., ¾" RG, 5-30°C, 0-30°dH
Water connection	Hot soft water: 1 pcs., 3/4" RG, 30-70°C, 0-3°dH
Water pressure	200-800 kPa
Outlet	Ø 110 mm P/S, wall or floor (to be specified)
Water consumption Urine bottle programme Urine bottle programme Extra-wash programme	2 wash phases, 16 I 3 wash phases, 24 I 5 wash phases, 40 I
Soap pump	Standard
Built-in detergent compartment	Room for 2 pcs. 5 I container
Pump capacity	130 l/min.
Dimensions wash chamber H x W x D	500 x 485 x 370 mm
Noise level	approx. 60 dB(A) (sound pressure)
Sealing class	Standard:IP 21 (not scavenging-proof)







iso **1400**1